



Septalloy is ideal for:

- Any posterior restoration
- Core build-up for crown preparation
- Core build-up with pre-fabricated posts and pins

Septalloy — the zinc free amalgam restorative that's a viable and cost effective alternative to composites.


Septalloy is the single-composition, non-gamma II, dental amalgam with a track-record for high performance and durability. A ternary alloy with high copper single composition, Septalloy exhibits excellent compressive strength and maintains marginal integrity over time. Its reliable consistency and easy handling make Septalloy the perfect solution to the demands of large posterior restorations as well as core build-ups for crown preparation and pre-fabricated posts.



From the manufacturers of Septocaine®

Easy-to-use — Easy-to-handle

Septalloy's convenient self-activating capsule provides one-step trituration without the use of accessories, to form a consistent, reliable amalgam alloy mix with superior handling characteristics. Color-coding of the reclosable capsules also enables easy differentiation between formulations and spill sizes.

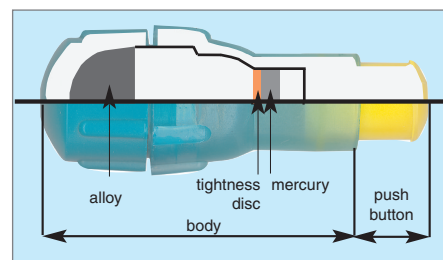

High Standards — High Quality — Exceptional Value


Septalloy is manufactured to the highest standards under the most strict quality control guidelines to deliver a superior alloy ideally formulated for any posterior restorations. That combined with a cost lower than most amalgams, makes Septalloy a great value.

Ternary Alloy with High Copper Single Composition

- Exhibits high compressive strength
- Resists marginal breakdown and clinical failure

Zinc Free

- Ensures superb handling
- Exhibits greater corrosion resistance

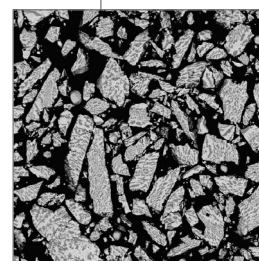

Physical Properties:

Alloy Composition

Ag:	46%
Sn:	31.3%
Cu:	22.7%

Hg/Powder Ratio:	1:2	Working Time:	3.5 minutes
Mercury Content:	54.5%	Carving Time:	8-10 minutes
Particle Shape:	Filings	Creep:	<1%
Mean Particle Size:	<40µm	Compressive Strength:	>350 MPa 24h:>
Recommended trituration time:	6-8 seconds	Dimension Variation:	-0.1 to +0.2%

- ▲ Tight Self-Activating Capsule
- ▲ Filing Shaped Particles
- ▲ Large Posterior Restoration

Ordering Information

ITEM NUMBER	DESCRIPTION
01-C1500 SEPTALLOY NG50	50 CAPSULES, 1 SPILL, CAP YELLOW, BODY GREEN
01-C1525 SEPTALLOY BULK NG50	500 CAPSULES, 1 SPILL, CAP YELLOW, BODY GREEN
01-C1550 SEPTALLOY NG50	50 CAPSULES, 2 SPILL, CAP GREEN, BODY GREEN
01-C1575 SEPTALLOY BULK NG50	500 CAPSULES, 2 SPILL, CAP GREEN, BODY GREEN



Visit www.septodontusa.com for more information.
800-872-8305 | septodontusa.com

