

The no-drip HEMA/Glutaraldehyde desensitizer in Glu/Sense™ goes where you want it to go.

Glu/Sense's unique SofNeedle™ tip and gel formula allows for precise placement. Guaranteed.

Glu/Sense is the only HEMA/Glutaraldehyde desensitizer on the market with a no-drip formula that provides control and accuracy.

Glu/Sense is a simple, one-step chairside procedure that treats and prevents dentinal hypersensitivity. With its 35% HEMA/5% Glutaraldehyde gel formulation, Glu/Sense acts in seconds to seal dentinal tubules, preventing the fluid shift, recognized as the primary cause of localized hypersensitivity.

- Unique gel formulation in a syringe, provides control and accuracy with no spills or drips on sensitive tissue
- SofNeedle™ Applicator
Foam-tip needle allows for faster, easier placement, without the burning of soft tissue
- Promotes adhesion
Minimizes microleakage around crowns and restorations
- Works with all bonding agents

Ordering is easy. Call 800-235-5862 or contact your preferred dealer. Visit centrixdental.com for more information.




CENTRIX
Glu/Sense™
HEMA/Glutaraldehyde Desensitizer



centrix®
Making Dentistry Easier.™

APPLICATIONS

- Under inlays, onlays, crowns or veneers
- Under amalgam fillings
- With bonding systems
- Where there are exposed dentinal tubules

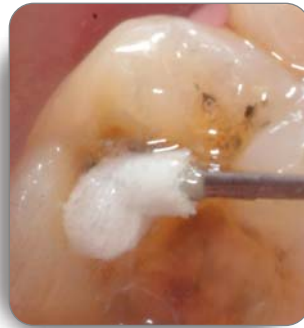

SofNeedle™ Applicator allows for precise placement of the gel formula.

PRODUCT DETAILS

- Rub on sensitive area for 30 seconds & allow to dry. Tooth does not need to be completely dry before application. No etching or light curing necessary. Avoid contact with soft tissue.
- Effect lasts 6 months or longer depending on patient. Permanent under restorations.
- Does not interfere with cements or bonding agents
- 2-year shelf life
- ADA CDT CODE 09910

Stop sensitivity before it starts.



ORDERING INFORMATION:

GLU/SENSE
REF 310122 6 x 1.0ml syringes
 + 60 SofNeedle Tips

SUGGESTED ACCESSORY:

REF 290038 SofNeedle 22ga tips, 144



Ordering is easy.

Call 1-800-235-5862
or contact your preferred dealer.

centrixdental.com



Our Guarantee — If you are not 100% satisfied, simply return the product for a complete refund. It's that easy.



centrix®

Making Dentistry Easier.™