

## **TIVA8 Dental**

HIGH THERMAL DISINFECTION WASHERS for Dental Instruments



Quality Products • Professional Excellence • Reliable Services • Time Saving Solutions

# Your Dental infection control success is backed by more than a century of global leadership experience

Tuttnauer's sterilization and infection control solutions are trusted by Dental and Medical clinics, laboratories, and hospitals in more than 140 countries. For over 95 years Tuttnauer has been passionate about delivering a complete sterilization solution for healthcare and life science industries.

We are proud to have installed over 350,000 sterilization and infection control products around the world.

Welcome to the TIVA8 Washer-Disinfector Designed to simplify your work and ensure the safety of your patients and your staff.





## **TIVA8 Dental 165 Liter Chamber**

High-end washer and disinfectant for washing, rinsing, disinfecting, and drying your dental instruments prior to sterilization.

UP TO 4 washing levels of total 8 DIN baskets capacity,

#### Improved Productivity and Efficiency



## **Intuitive User Experience**

- A full control and all-in-one operation from a single touch panel with color display
- Graphical report of washing phases including time and temperature
- Double-layered tempered glass door for visual monitoring
- Spray Arm monitoring for continuous washing validation (optional)
- Multi color chamber LED light (White, Green, Red) operates during loading, unloading and running the program (optional)



#### **Eco-Friendly**

- Cycle time reduction thanks to a preheated water tank (only for HS model)
- 2 dosing pumps with flowmeter for precise consumption of acid and alkaline detergents
- Integrated water softener for reducing decontamination
- Built-in steam condenser for excess water and steam



#### **Enhanced Traceability**

- Comply with regulation by saving washing cycles data
- USB and RS232 connectivity
- Printer for data documentation (optional)



#### **Premium Performance**

- Instruments are completely dried thanks to powerful electric blowers and heaters with a HEPA H14 filter
- Up to 3 water connections for flexible operation and for optimal washing needs (cold, hot, demineralized)
- Up to 40 washing programs
- AISI 316L Stainless Steel

# TIVA 8 UNDERCOUNTER MANUAL DOOR

#### UP TO 4 washing levels of total 8 DIN baskets capacity, 165 Liter chamber

The TIVA 8 washer-disinfector is equipped with a manual double glass door and is best suited for clinics and dental practices. It can be loaded with up to 4 washing levels with a large variety of washing trolleys and accessories. Available trolleys for DIN baskets, containers, mini-invasive equipment, anesthesia kits, bottles, shoes and dental handpieces.

#### TIVA8-1M

Narrow Design 165 Liter chamber Dimensions (WxDxH mm) 600x650x860



#### TIVA8-WD-1M / TIVA8-HS-1M

Wide design 165 Liter Chamber Dimensions (WxDxH mm) 900x650x860





Selectable configurations:

Detergent storage cabinet /Pre-heating tank for reduction of cycle time + Detergent cabinet and demineralized water connection

| Model                            | Wash Chamber<br>dimension<br>WxHxD (mm) | Chamber<br>Vol.<br>(Liters) | Capacity                     | Overall<br>Dimensions<br>WxHxD (mm) | Weight<br>(kg)   | Power Supply   | Power<br>Consumption | Item supplied                        |
|----------------------------------|---|-----------------------------|------------------------------|-------------------------------------|------------------|--|----------------------|--------------------------------------|
| TIVA 8 -1M                       |   |                             | Up to 4<br>washing<br>levels | 600x860X650                         | 100 kg           | 230V, 1 phase, 50Hz<br>200V, 1 phase, 50/60Hz<br>380/415V, 3 phase | 6.3 kW               | Drain tube,<br>Water supply<br>pipes |
| TIVA 8-WD-<br>1M<br>TIVA 8-HS-1M | 550x600x500                             | 550x600x500 165L (8 DIN     | (8 DIN<br>baskets)           | 900x860X650                         | 115 kg<br>120 kg | +Neutral, 50Hz<br>220V, 3 phase, 60Hz<br>50/60Hz<br>50/60Hz        | 7 kW<br>8.8 kW       | Drain tube,<br>Water supply<br>pipes |

#### **Features**

- Manual double glass door
- Washing chamber made of high quality AISI 316L
- Touch panel with color graphic display
- A total of 40 programs: Preset and programmable procedures
- Air drying prefilter + HEPA H14 filter
- Preheated tank for cycle time reduction (only for HS model)
- 2 dosing pumps with flowmeter for acid and alkaline
- detergents (2 additional pumps optional)

- 3 heating elements
- USB/ Printer for data documentation
- Built-in steam condenser

#### **Optional**

- Conductivity measurements to check rinse water quality before disinfection phase
- Demi water boiler (only for HS version)
- Multi color LED chamber light: Operates during loading, unloading and running the program.
   The LED vary in color depends on machine status: White: open door | Green: cycle ends | RED: error
- Spray arm monitoring continuous validate rotation of spray arm to guarantee that the load is currently washed
- Stainless steel stand for ergonomic loading and unloading (not installed under-counter)
- Water connection: Cold water softener, purified water (TIVA 8-1M, TIVA 8-WD-1M) Cold water softener (TIVA 8-HS-1M)
- Thermal Printer Graphic report of washing phases including time and temperature


#### **Drawings:**

#### **TIVA 8-1M**

#### TIVA 8-WD-1M / TIVA 8-HS-1M

60 cm narrow design

90 cm wide design

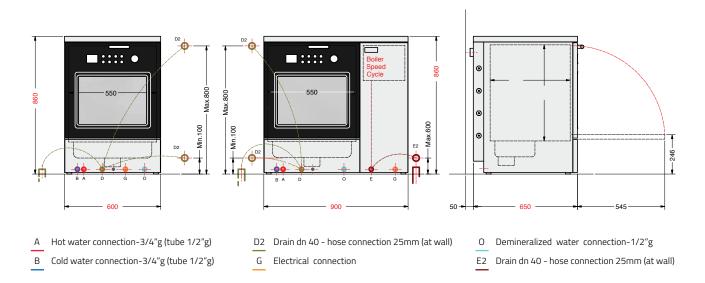

## **SAVIT**

#### RACKS

| Code      | Image  | Description   |                |  |
|-----------|--|---|----------------|--|
| 119905148 |  | TH148 LOWER WASH RACK   |                |  |
|           |  | Washing rack positioned on the lower level, basic version For various inserts and accessories.  |                |  |
|           |  | Internal loading dimensions with upper rack WxDxH   | 500x475x230 mm |  |
|           |  | Internal loading dimensions without upper rack WxDxH  | 500x475x470 mm |  |
|           | 0  | Rack external dimensions WxDxH  | 525x502x48 mm  |  |
|           |  | Empty rack weight   | 2,5 Kg         |  |
|           | Configuration example with:<br>2x A085 + 1x A053 |   |                |  |
| 119905156 |  | TH156 LOWER WASH RACK 2 CONN.  Washing rack positioned on the lower level, with 2 connection Left-hand side free for various inserts (475x215 mm).  Active drying of turbines, hand-pieces and contraangles. For various inserts and accessories. | 1              |  |
|           |  | Internal loading dimensions with upper rack WxDxH   | 500x430x230 mm |  |
|           |  | Internal loading dimensions with upper rack and A076 WxDxH  | 360x460x230 mm |  |
|           |  | Internal loading dimensions without upper rack WxDxH  | 500x430x470 mm |  |
|           |  | Rack external dimensions WxDxH  | 525x525x72 mm  |  |
|           |  | Empty rack weight   | 3 Kg           |  |
|           | Configuration example with:                      |   |                |  |
|           | 1x A085 + 1x A053 + 1x A076 + 2x A075            |   |                |  |

#### Code Image

#### **Description**

119905157



#### TH157 LOWER WASH RACK 4 CONN.

Washing rack positioned on the lower level, with 4 connection

Active drying of turbines, hand-pieces and contra-angles. For various inserts and accessories.

| Internal loading dimensions with upper rack WxDxH          | 500x430x230 mm |
|--|----------------|
| Internal loading dimensions with upper rack and A076 WxDxH | 220x460x230 mm |
| Internal loading dimensions without upper rack WxDxH       | 500x430x470 mm |
| Rack external dimensions WxDxH                             | 525x525x72 mm  |
| Empty rack weight  | 3,5 Kg         |



Configuration example with: 2x A076 + 4x A075 + 1x A085


#### TH151 UPPER WASH RACK W/S.ARM


Washing rack positioned on the upper level, with spray arm

For various inserts and accessories.

| Internal loading dimensions WxDxH | 480x460x185 mm |
|-----------------------------------|----------------|
| Rack external dimensions WxDxH    | 532x524x110 mm |
| Empty rack weight                 | 3,5 Kg         |



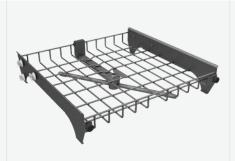

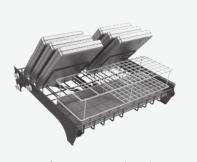


Configuration example with: 2x A065 + 1x A053

## **TIVA8**

#### **RACKS**

## Code **Image Description** 119905149 TH149 UPPER WASH RACK 2 CONN. W/S.ARM Washing rack positioned on the upper level, with 2 connection and spray arm Left-hand side free for various inserts (475x215 mm). Active drying of turbines, handpieces and contraangles. For various inserts and accessories. Internal loading dimensions with A076 WxDxH 350x460x185 mm Internal loading dimensions without A076 WxDxH 480x430x185 mm Rack external dimensions WxDxH 532x524x110 mm Empty rack weight Configuration example with: 1x A053 + 1x A063 + 1x A076 + 2x A075 119905153 **TH153 UPPER WASH RACK 4 CONN.** Washing rack positioned on the upper level, with 4 connection and spray arm Active drying of turbines, hand-pieces and contra-angles. For various inserts and accessories. Internal loading dimensions with 2x A076 WxDxH 220x460x185 mm Internal loading dimensions without A076 WxDxH 480x430x185 mm Rack external dimensions WxDxH 532x524x110 mm 5 Kg Empty rack weight Configuration example with: 2x A076 + 4x A075 + 1x A085

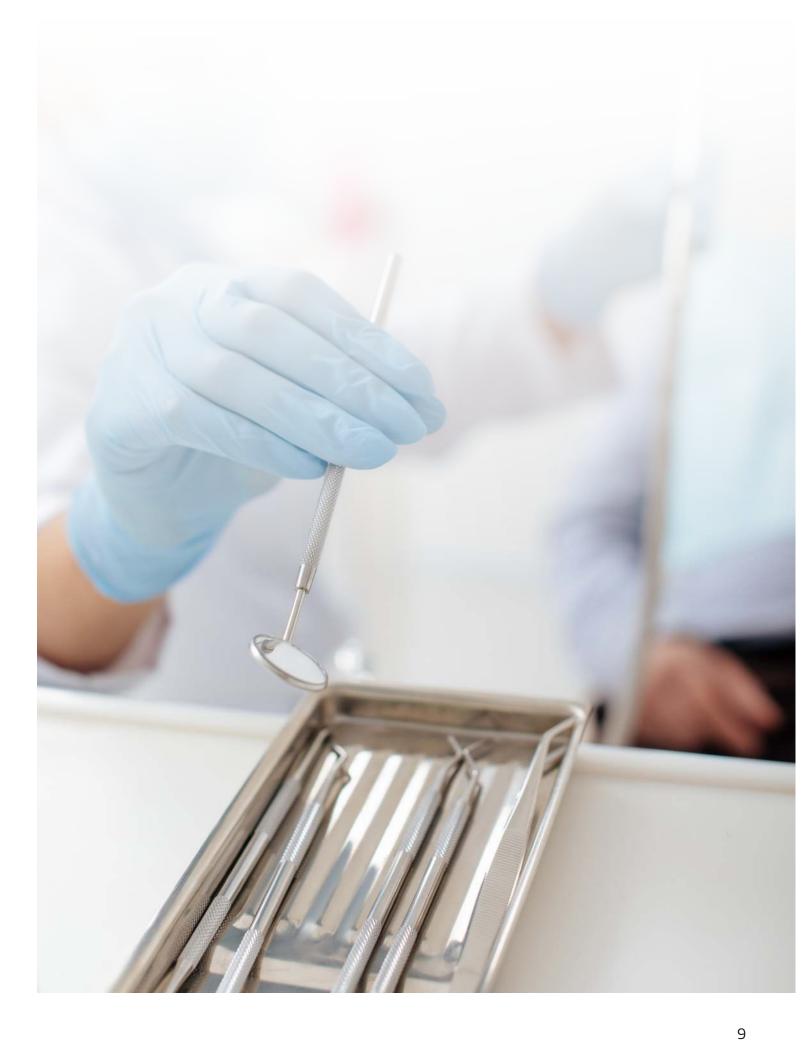

## **SAVIT**

#### ACCESSORIES

| Code        | Image | Description  |
|-------------|-------|--|
| P00001039   |       | A076 INSERT FOR HANDPIECES HOLDER Insert for 2 holders Suitable for use on all racks To be used in combination with:  • A084 • A075 Insert external dimension WxDxH 104x455x77 mm  |
| P00001048   |       | A075 HOLDER 8 POS.  Preassembled holder with 8 positions for handpiece holder Includes: S/S FILTER  To be used in combination with:  • A076 Possibility of mounting up:  • 8x A066  • 8x A079  • 8x A127  • 8x A135                                      |
| P00001098   |       | A084 HOLDER 12 POS.  Preassembled holder with 12 positions for handpiece holder Includes: S/S FILTER  To be used in combination with:  • A076  Possibility of mounting up:  • 12x A066  • 12x A127  • 6x A079 + 6x A066  • 6x A067 + 6x A135  • 12x A135 |
| ATS238-0412 |       | A079 SUPPORT HANDPIECES HOLDER To be used in combination with:  • A075 • A084  |
| ATS238-0411 |       | A077 HANDPIECES ADAPTER Ø19-21  Adapter for handpieces with dimension from Ø 19 to Ø 21 mm To be used in combination with:  • A079   |
| ATS238-0410 |       | A078 HANDPIECES ADAPTER Ø14-17 Adapter for handpieces with dimension from Ø 14 to Ø 17 mm To be used in combination with:  • A079  |
| ATS238-0387 |       | A067 HANDPIECES CAP To be used in combination with:  • A077  • A078  |
| ATS238-0438 |       | A080 SUPPORT FOR 8 SCALER TIPS  To be used in combination with:  ■ A077  |
| D00564000   |       | A139 SUPPORT FOR SILICONE CONN.  To be used in combination with:  • A075  • A084  • A138   |

| Code        | Image  | Description  |
|-------------|--|--|
| D00291500   | 6  | A138 SILICONE CONN.  To be used in combination with:  A139   |
| D00109200   | TEE  | A140 CAP FOR SILICONE CONN.  To be used in combination with:  A138   |
| ATS238-0386 |  | A066 NOZZLE FOR ASPIRATOR TIPS  To be used in combination with: ■ A075 ■ A084  Nozzle external dimension Ø 4 x H80 mm                        |
| D00606500   |  | A142 ADAP. HOSE 1/4" TO G1/8" M. THREAD Stainless steel hose male holder adapter DN 5 and DN 6   |
| 467010016   |  | A127 ADAP. MALE LUER TO G1/8" M. THREAD Stainless steel male Luer-Lock adapter   |
| 467010018   | 9  | A129 PVDF FEMALE LUER TO 3/16" HOSE HOLDER PVDF hose holder DN 5 to female Luer-Lock   |
| 467010019   |  | A131 PVDF MALE LUER TO 3/16" HOSE HOLDER PVDF hose holder DN 5 to male Luer-Lock   |
| D00341100   |  | A149 MALE LUER TO 1/4" HOSE HOLDER Stainless steel hose holder DN 6 to male Luer-Lock  |
| 467010020   | 9  | A132 PVDF FEMALE LUER TO 1/4" HOSE HOLDER PVDF hose holder DN 6 to female Luer-Lock  |
| 467010021   |  | A133 PVDF MALE LUER TO 1/4" HOSE HOLDER PVDF hose holder DN 6 to male Luer-Lock  |
| 467010022   | P  | A134 PVDF ELBOW FITTING 1/4" HOSE HOLDER PVDF elbow hose holder DN 6   |
| 467010023   |  | A147 PVDF FEMALE LUER CAP PVDF female Luer-Lock cap  |
| 467010024   |  | A148 PVDF M. LUER CAP PVDF male Luer-Lock cap  |
| D00304700   | The state of the s | A144 FITTING 3/16" HOSE HOLDER Stainless steel pass trought hose holder DN 5   |
| ATS238-0439 |  | A081 LUER-LOCK CONNECTOR KIT  Suitable for use with racks:  TH121 TH123  Specifications: Female Luer-Lock connection at the end of the pipe. |

## **SAVIT**

#### ACCESSORIES

| Code        | Image | Description   |
|-------------|-------|---|
| ATS238-0434 |       | A135 CAP SCREW G1/8"  To be used in combination with:  • A075  • A084   |
| ATS238-0385 |       | A069 NEMA HANDPIECES ADAPTER  It allows to have the perfect inside channels cleaning through water circulation. To be used in combination with:  • A077  • A078   |
| ATS238-0384 |       | A073 BASKET WITH LID  Rinsing basket with fixed lid for small accessories Made of high-quality stainless-steel wire for storing small parts and accessories during the cleaning process.  Basket external dimension WxDxH 80x80x40 mm |
| ATS238-0381 |       | A068 INSERT FOR 10 DENTAL IMPRESSION Suitable for use on all racks Insert external dimension WxDxH 184x131x130 mm   |
| ATS238-0382 |       | A071 INSERT FOR 7 EXTRACTION PLIERS Suitable for use on all racks Insert external dimension WxDxH 184x131x122 mm  |
| ATS238-0306 |       | A055 INSERT FOR 9 TRAYS Suitable for use on all racks Insert external dimension WxDxH 455x245x127 mm  |
| ATS238-0376 |       | A063 INSERT 8 CASSETTES Insert with 8 positions for trays or cassettes with width max 35 mm Suitable for use on all racks Insert external dimension WxDxH 180x370x105 mm  |
| ATS238-0380 |       | A070 INSERT 7 POS. Insert with 7 positions for trays or cassettes with width max 40 mm Suitable for use on all racks Insert external dimension WxDxH 180x370x105 mm   |
| ATS238-0444 |       | A088 INSERT 3 POS. Insert with 3 positions for trays or cassettes with width max 40 mm Suitable for use on all racks Insert external dimension WxDxH 180x360x135 mm   |
| ATS238-0379 |       | A065 INSERT 3 CASSETTES Insert with 3 positions for trays or cassettes with width max 35 mm Suitable for use on all racks. Insert external dimension WxDxH 170x255x140 mm   |

| Code        | Image                            | Description   |
|-------------|----------------------------------|---|
| ATS238-0447 |                                  | A087 BASKET SPACERS  To be used in combination with:  • A082  • A083  • A085  Specifications: used to stock 2 baskets A083 or 2 baskets A082 or 2 baskets A085.   |
| D00552800   |                                  | A053 INSERT FOR VERTICAL INSTRUMENTS Suitable for use on all racks. Insert external dimension WxDxH 445x165x120 mm Space requirements into the rack: 1/3  |
| 000552600   | Space requirements into the rack |   |
| D00594900   |                                  | A085 1/3 BASKET  Dental basket 1/3 For various inserts and accessories. Suitable for use on all racks.  Insert external dimension WxDxH 455x165x40 mm  Space requirements into the rack: 1/3  |
| 000394900   | Space requirements into the rack |   |
| ATS238-0305 |                                  | A054 BASKET WITH 2 SILICONE SUPPORT Stainless steel basket with two silicone supports. Suitable for use on all racks. Basket external dimension WxDxH 200x280x30 mm   |
| ATS238-0103 |                                  | A136 STAINLESS STEEL STAND TIVA8 Stand for TIVA8-1M (and variants) provide an ergonomics height of 146 cm, when the machine is not installed undercounter. Suitable for use as a detergent storage compartment Stand external dimension WxDxH: 600x585x600 mm.      |
| ATS238-0104 |                                  | A137 STAINLESS STEEL STAND L VERS. Stand for TIVA8-WD-1M (and variants) provide an ergonomics height of 146 cm, when the machine is not installed undercounter. Suitable for use as a detergent storage compartment Stand external dimension WxDxH: 900x585x600 mm. |
| P00001065   | Listersess                       | KIT DISPENSER VALIDATION Kit for calibrating dispensers.  |


## SPARE PARTS



## **TIVA8**

#### LAYER CONFIGURATION SHCEME


#### **Useful space dimensions**

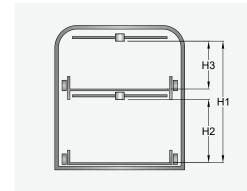

## DOUBLE LEVEL WITHOUT HOLDER AND SUPPORT FOR HANPIECES ACCESSORIES

Configuration example: TH151+ TH148


| H1     | H2     | H3     |  |
|--------|--------|--------|--|
| 470 mm | 230 mm | 185 mm |  |

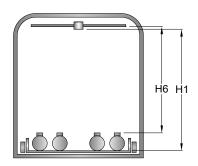

## LOWER LEVEL WITH HOLDER AND SUPPORT FOR HANPIECES ACCESSORIES

Configuration example: TH157

| H1     | Н6     |  |
|--------|--------|--|
| 470 mm | 415 mm |  |

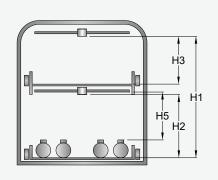


## LOWER LEVEL WITH HOLDER AND SUPPORT FOR HANPIECES ACCESSORIES

## UPPER LEVEL WITHOUT HOLDER AND SUPPORT FOR HANPIECES ACCESSORIES

Configuration example: TH157 + TH148

| H1     | H2     | НЗ     | H5     |
|--------|--------|--------|--------|
| 470 mm | 230 mm | 185 mm | 174 mm |

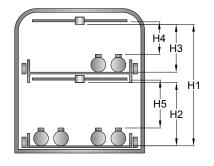

## DOUBLE LEVEL WITH HOLDER AND SUPPORT FOR HANPIECES ACCESSORIES

Configuration example: TH157 + TH149

| H1     | H2     | Н3     | H4     | H5     |
|--------|--------|--------|--------|--------|
| 470 mm | 230 mm | 185 mm | 130 mm | 174 mm |

| IOTES |  |
|-------|--|
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |

#### **Washing & Disinfection Indicators**



WTL198-0062



WTL198-0061



WTL198-0065

Cleaning Indicator Holder



WTL198-0063

## Your clinic's infection control success is backed by 95 years of experience and global leadership

Renowned global leaders in infection control solutions with more than

95 years experience

Products distributed to clinics, hospitals and research institutions in more than

countries

#### International Sales and Marketing

E-mail: info@tuttnauer.com www.tuttnauer.com

#### Tuttnauer Europe b.v.

Hoeksteen 11, 4815 PR PO Box 7191, 4800 GD Breda The Netherlands Tel: +31 765 423 510 Fax: +31 765 423 540

E-mail: info@tuttnauer.nl

#### Tuttnauer USA Co.

25 Power Drive, Hauppauge, NY 11788 Tel: +800 624 5836, +631 737 4850 Fax: +631 737 0720 E-mail: info@tuttnauerUSA.com

Distributed by:



in 5 A D







Learn from our experts Join our blog tuttnauer.com/blog